


**ACCESSORIES**

**Rain Cover**

**# 996625**

Item	Description	Qty
1	Rain Cover	1

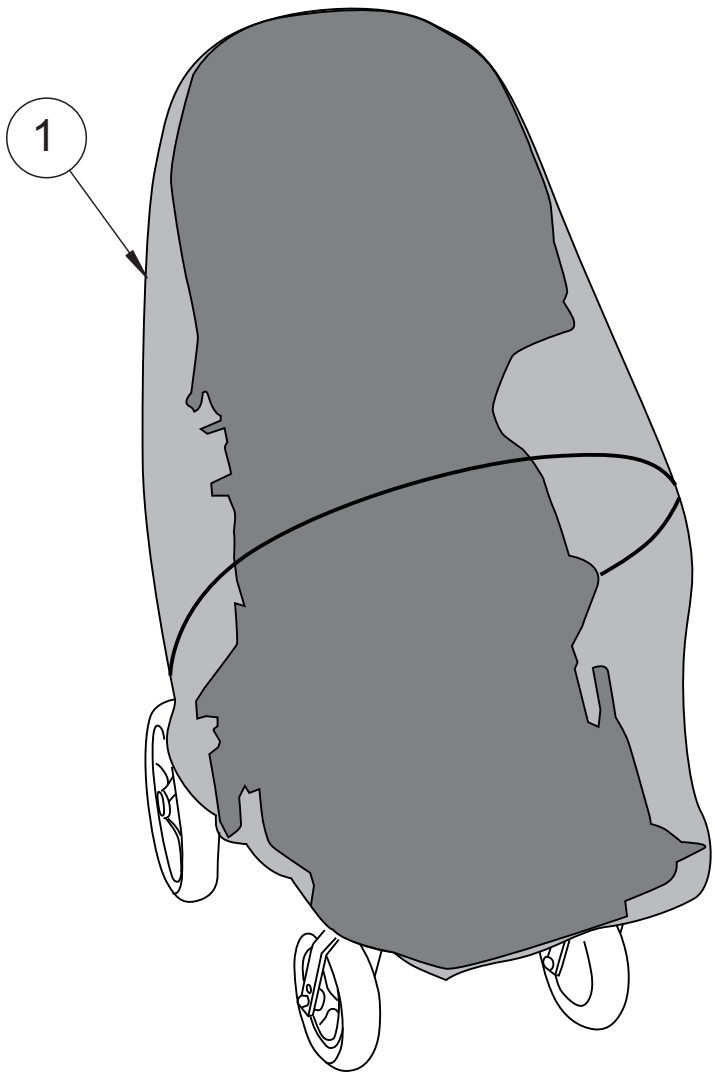

If you have any questions about attaching the Rain Cover, please call our Sales Department at 888 CONVAID (266-8243) or +1 310 618-0111.

## INSTRUCTIONS

---

### Rain Cover

1. Take the top of the rain cover and wrap it around the top of the canopy
2. Arrange the rain cover around the wheelchair and pull the bottom of the rain cover to wrap around the footplate
3. Find the two Velcro® straps on the right and left side of the rain cover.
4. Secure rain cover to wheelchair by wrapping Velcro® straps around each side of wheelchair frame.

**NOTE: For smaller size wheelchairs, the Velcro® straps may connect to each other behind the chair for a tighter fit.**

Available in two sizes:

Small – fits chair sizes 10-14

Large – fits chair sizes 16-18