

In Stock

Convaid

CARROT


3

Special Needs Child Restraint System

- Uniquely designed to grow from ages 3-15 (30 - 108 lbs. / 37" - 60")
- Modular extensions provide extensive growth and positioning where needed - add length to the seat, height to the back...or both!
- High density foam seat conforms to the natural contour of the body to provide maximum comfort & positioning
- "Free angle" recline opens the hip area to provide added comfort, and promotes head & trunk control
- Reclines flat for diaper changes
- Ergonomic head support provides optimal positioning and side impact protection
- Firm lateral supports maintain midline positioning of the trunk for increased postural control
- Depth adjustable pommel separates the thighs to improve hip alignment
- Positioning inserts contour the shape of the Carrot to give your child a custom fit
- Body contouring design protects more securely against side impacts
- LATCH equipped (Lower Anchors and Tethers for Children)
- Exceeds U.S. Federal Motor Vehicle Safety Standards FMVSS 213 & 302
- Certified for aircraft travel
- 5 year warranty - one of the best in the industry
- Sold in over 16 countries - top selling Restraint System in Europe!



Recline Adjustment
90° - 180°

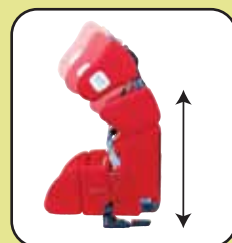

Back Height
11.8" - 22.8"



Seat Depth
9.8" - 22.3"

 **Convaid**

1-888-convaid • www.convaid.com

Features

Carrot 3 Special Needs Booster



Head Support

Ergonomically designed Head Support provides optimal support and protects the head from side impacts and centrifugal forces.



Shoulder Wings

Expandable Shoulder Wings protect the body from side impacts and centrifugal forces.



Firm Lateral Supports

Maintains midline positioning of the trunk for increased postural control.



Chest Support

The Chest Support holds the occupant in a stationary upright position.



Pommel/Abductor


Separates the thighs to improve hip alignment and stabilize sitting posture.



LATCH Connector System

LATCH connectors (top tether strap & lower anchors) securely install the Carrot into the vehicle providing the safest seating possible. Note: The Carrot 3 must be used with a vehicle lap and shoulder belt. Integrated 5-point harness is for positioning only.

Growth Adjustment Chart

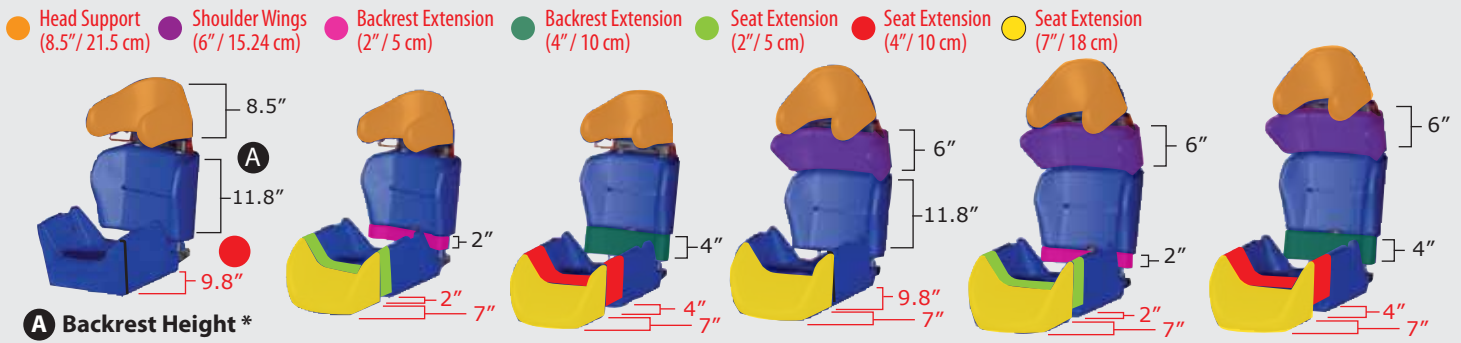

in (A)	11.8	13.8	15.8	17.8	19.8	21.8
cm (A)	30	35	40	45.2	50.3	55.4

in (B)	9.8	11.8	13.8	9.8	11.8	13.8
cm (B)	24.9	30	35	24.9	30	35

* Backrest Height can be raised an additional 1" in addition to the 4" seat extension. ^ Seat Depth can be lengthened an additional 1.5" without using a rodded extension.

Optional Accessories

Shoulder Wings	Rodded Seat Extension - Large	Rodded Seat Extension - Short	Support Tray (Crash Tested)	Footrest (S/M/L)
Lateral Support	Seat Extension (2" / 4")	Backrest Extension (2" / 4")	Pommel	Top Tether Strap (Required)