


## Rodeo Tilt Features

TILT

- 10° - 45° of Adjustable Tilt  
*Provides pressure relief and upper body support for individuals with insufficient trunk stability*
- 90° - 105° Seat-to-Back Angle Adjustment  
*Provide better positioning for individuals with limited ability to bend at the waist*
- Firm Lateral Thoracic Supports (option)  
*Maintain midline positioning of the trunk for better respiration and digestion*
- In-Depth Positioning Options  
*Provide maximum comfort and optimal positioning*
- Size Adjustable Frame  
*Extends the life of the chair for a growing child*
- Lightweight Frame  
*Is easy to lift and carry - one of the lightest tilt-in-space wheelchairs available!*
- Tilt-in-Space Base with Unique Compact Fold  
*Folds compactly for easy transport and storage - fold and go!*
- Rigorous Product Testing  
*Ensures maximum safety, strength, stability, and endurance*
- Crash Tested Transit Models  
*Adhere to national safety standards, and allow safe transport in bus or van*

RODEO

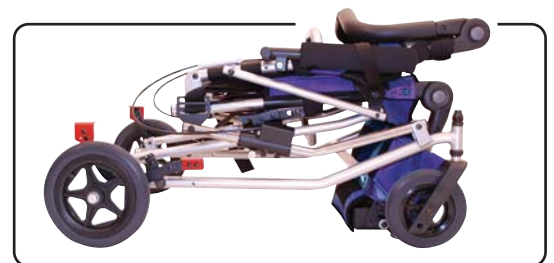

compact folded  
**RODEO**

# Rodeo Tilt Specifications

Measurement	Chair Model	Seat Width*	Seat Depth*	Seat Back Height	Seat to Footrest*	Footplate Adjustment Angle	Seat to Floor	Seat-to-Back Angle	Tilt Adjustment (from vertical)	Headrest Extension	Shoulder Strap Height	Weight Capacity Regular/Transit	Weight of Chair <sup>Δ</sup>
in.	<b>RD 12</b>	12	11 - 14	23	8 - 13	15°	20	90°-105°	10° - 45°	6	19 - 22	75 lbs / 66 lbs	28.3 lbs
cm	<b>RD 12</b>	30.5	28 - 35.5	58.4	20 - 33		51			20	48 - 56	34 kg / 30 kg	10.5 kg
in.	<b>RD 14</b>	14	13 - 16	25	9.75 - 15	15°	21	90°-105°	10° - 45°	8	18.5 - 24	100 lbs / 100 lbs	29.5 lbs
cm	<b>RD 14</b>	35.5	33 - 41	63.5	25 - 38		53			20	47 - 61	45.5 kg / 45.5 kg	12 kg
in.	<b>RD 16</b>	16	15 - 18	27	12 - 17.5	15°	23	90°-105°	10° - 45°	8	19 - 25.5	170 lbs / 170 lbs	35.6 lbs
cm	<b>RD 16</b>	40.5	38 - 46	69	30 - 44		58			20	48 - 64.8	77 kg / 77 kg	14 kg

in. - Inches, cm - Centimeters \* Plus or minus .5 inch (1.27 cm) due to fabric. Δ - Weight of chair is without front rigging.

Use only Convoid options and parts on Convoid products. Convoid parts are not interchangeable with other manufacturer's products. Replace any worn parts immediately.

Note: Convoid maintains a policy of continual product improvement; therefore, we reserve the right to make changes without notice.

Measurement	Chair Model	Overall Width	Overall Height	Overall Length	Folded Length	Folded Width	Folded Height
in.	<b>RD 12</b>	20	41	31	33	20	17
cm	<b>RD 12</b>	51	104	79	84	51	43
in.	<b>RD 14</b>	22	44	34	34	22	17
cm	<b>RD 14</b>	56	112	86.4	86.4	56	43
in.	<b>RD 16</b>	24.5	45	39	34.5	24.5	19
cm	<b>RD 16</b>	62	114	99	87.5	62	48

## Limited Lifetime Warranty

Frame and X-Braces.....Lifetime of original retail buyer  
 Fabric and webbing.....One year  
 Other components.....One year

This warranty does not cover normal wear and tear, or damage caused by accident or misuse.

## Upholstery Color Options

### Upholstery Colors Removable, Washable Fabric



One Color Rodeo

Two Color Rodeo

U.S. Patent 6,345,835

This product may be covered by the following U.S. patent - 6,099,075  
 Additional patents pending.

