

## Safari Tilt Features


# TILT

- 5° - 45° of Adjustable Tilt  
*Provides pressure relief and upper body support for individuals with insufficient trunk stability*
- 90° - 135° of Adjustable Recline  
*Opens the seat-to-back angle for better digestion, improved respiration, and pressure relief*
- Support Cushion  
*Promotes stability and posture control, and waterfall front eliminates pressure points at the back of the knees*
- Elevating Legrests (option)  
*Allows for changes in position and promotes better circulation*
- Tilt-in-Space Wheelchair with Unique "Umbrella Style" Fold  
*Folds compactly for easy transport and storage – fold and go!*
- In-Depth Positioning Options  
*Provide maximum comfort and optimal positioning*
- Size Adjustable Frame  
*Extends the life of the chair for a growing child*
- New Wider Wheels  
*Provide added stability and a comfortable ride on uneven terrain*
- Crash Tested Transit Models  
*Adhere to national safety standards, and allow safe transport in bus or van*

# SAFARI



One-Piece Footplate (Standard)

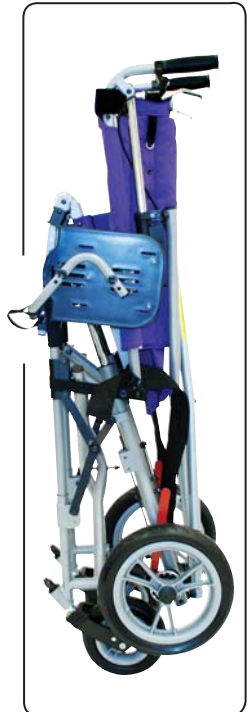

compact folded  
**SAFARI TILT**

**SFT12/14/16:** 7.5" x 2" Front, 11.5" x 2.5" Quick Release Rear Wheels  
**SFT18:** 8" Front, 11" Rear

Positioning options and accessories sold separately.

# Safari Tilt Specifications

Measurement	Chair Model	Seat Width*	Seat Depth	Seat Back Height	Seat to Footplate*	Seat to Floor*	Armrest to Seat	Seat to Back	Recline Adjustment	Tilt Adjustment (from vertical)	Headrest Extension	Shoulder Strap Height
in.	<b>SF 12</b>	12	10 - 14	23	11 - 13	23.5	8 - 11	90° - 95°	3°, 22°, 45°	5° - 45°	5	8.5 - 20
cm	<b>SF 12</b>	30.5	26 - 36	58	28 - 33	59.7	21 - 27				13	22 - 51
in.	<b>SF 14</b>	14	12 - 16	22	11 - 13	26	8 - 10.5	90° - 95°	1°, 19°, 40°	5° - 45°	8	9.5 - 19.5
cm	<b>SF 14</b>	35.5	30 - 41	56	28 - 33	66	20 - 26				20	24 - 50
in.	<b>SF 16</b>	16	14 - 18	22	12.5 - 14.5	26	8 - 11	90° - 95°	1°, 19°, 40°	5° - 45°	8	10 - 20
cm	<b>SF 16</b>	40.5	36 - 46	56	31.8 - 36	66	21 - 27				20	25 - 51
in.	<b>SF 18</b>	18	16 - 20	21.5	10 - 16	23	7.5 - 10	90° - 95°	-----	5° - 45°	8	12.5 - 18.5
cm	<b>SF 18</b>	35.5	40.5 - 51	55	25 - 41	58	19 - 25				20	32 - 47

in. - Inches, cm - Centimeters \* Plus or minus .5 inch (1.27 cm) due to fabric.

Measurement	Chair Model	Overall Width	Overall Height	Overall Length	Folded Length	Folded Width	Folded Height	Weight Capacity Standard/Transit	Weight of Chair $\Delta$ *
in.	<b>SF 12</b>	21.5	44.5	36.5	48	16.5	14.5	75 lbs / 66 lbs	34 lbs
cm	<b>SF 12</b>	54.6	113	92.7	122	42	36.8	34 kg / 30 kg	15.5 kg
in.	<b>SF 14</b>	24.3	45.5	39	50	14.8	14	100 lbs / 100 lbs	37 lbs
cm	<b>SF 14</b>	61.6	115.6	99	127	37.5	35.6	45.5 kg / 45.5 kg	16.8 kg
in.	<b>SF 16</b>	25.5	45.5	42.5	51	16.25	14	170 lbs / 170 lbs	41 lbs
cm	<b>SF 16</b>	64.8	115.6	108	129.5	41.3	35.6	77 kg / 77 kg	18.6 kg
in.	<b>SF 18</b>	26	44.5	41.5	51.5	15	14	250 lbs / 170 lbs	45 lbs
cm	<b>SF 18</b>	66	113	105.5	131	38	36	114 kg / 77 kg	20.5 kg

$\Delta$  - Weight of chair is without front rigging. \* Standard model chair with recline option add 3 lbs. (1.4 kg) to weight of chair.

Use only Convoid options and parts on Convoid products. Convoid parts are not interchangeable with other manufacturer's products. Replace any worn parts immediately.

Note: Convoid maintains a policy of continual product improvement; therefore, we reserve the right to make changes without notice.

## Upholstery Colors



## Limited Lifetime Warranty

Frame and X-Braces.....Lifetime of original retail buyer  
 Fabric and webbing.....One year  
 Other components.....One year

This warranty does not cover normal wear and tear, or damage caused by accident or misuse.



U.S. Patent 6,206,405  
 This product may be covered by the following U.S. patent - 6,099,075 and Japan 3072085  
 Additional patents pending.

Rev 8-4-10

## Seating Options



### Support Cushion

Flat Solid Cushion (Standard)

Promotes stability and posture control. Waterfall front eliminates pressure points at back of knees.

E2601 - General use wheelchair cushion.



### Position Cushion

Lateral Pelvis, Lateral and Medial Thigh Support/Anti-Thrust Cushion

Medial thigh support positions legs out of adduction and into better alignment.

E2607 - Skin protection and positioning wheelchair cushion.



### Align Cushion

Lateral Pelvis/Thigh Support Cushion

Contoured edges provide lateral pelvic and thigh positioning.

E2607 - Skin protection and positioning wheelchair cushion.